

Acquipix User Manual

Version 0.5, 20210805

By Dr. Jorrit Montijn

Cortical Structure and Function group

Netherlands Institute for Neuroscience

E-mail: j.montijn@nin.knaw.nl

Contents

Overview.....	3
Installation instructions.....	4
Before you start.....	4
Installing the external libraries for the recording computer.....	4
Installing the external libraries for preprocessing your data	4
Installing Acquipix and its dependencies	4
Workflow summary	5
User guide.....	6
Running an experiment	6
Using the online analysis tools	7
Preprocessing your data.....	8
Setting the variables and parameters	8
Custom data-preprocessors of NI channels	10
Compiling the library and pre-processing the data.....	11
Adjusting the probe coordinates and exporting the data	12
Exporting data to the Acquipix (AP) structure format	13
Exporting data to the Neurodata Without Borders (NWB) format.....	13
Technical descriptions.....	14
The Acquipix (sAP) data structure format.....	14
Overview.....	14
The cellBlock structure	14
The sStimObject substructure	15
The sCluster structure	16
Experiment scripts.....	17
Overview.....	17
Description of inputs and structure of an experiment script.....	17
Structure of the experiment script & output file	19
Integrating your experiment file into Acquipix	21
Online analysis functions.....	21
Overview.....	21
Troubleshooting	22

Overview

Multi-stream data recorded while performing a visual or optogenetic experiment, can be difficult and laborious to process. If your data includes stimulation, multi-channel spiking, LFPs, running speed, and pupil tracking, you will likely already have to deal with 5 independent data streams that need to be synchronized (fig. 1). Moreover, keeping track of all these different files can be a pain. That's why we created the Acquipix repository, which does all the synchronizing and synthesizing for you. All our code is modular, open source and fully customizable, but we made sure you can also run everything out-of-the-box using only graphical user interfaces without having to write a single line of code.

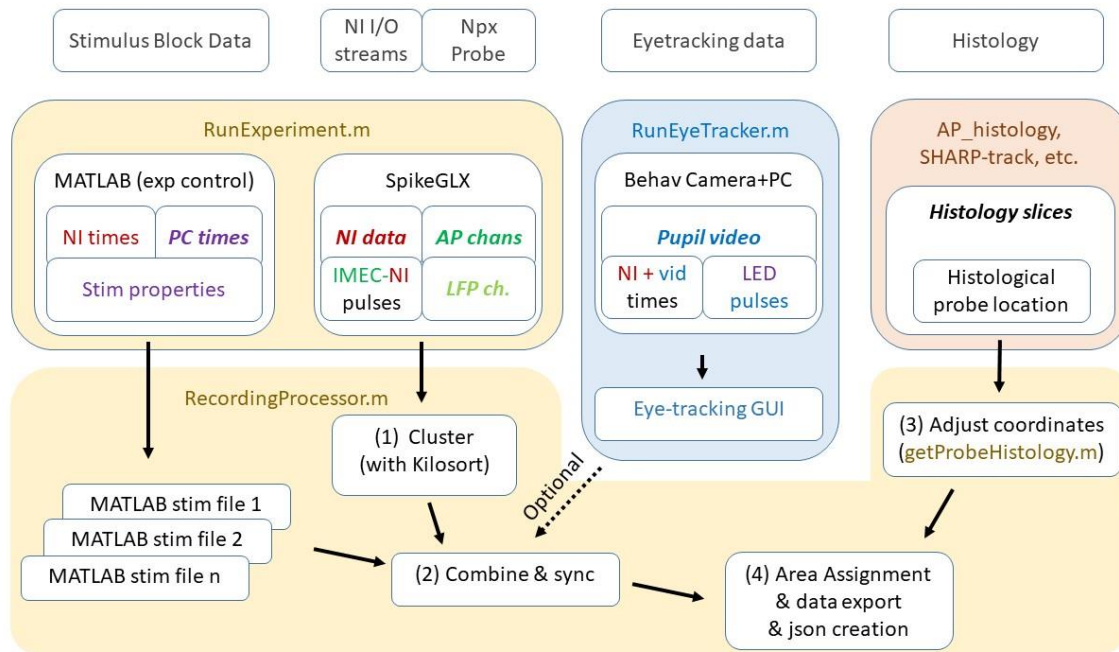


Figure 1. Technical data flowchart showing a single recording. First column from left (*MATLAB (exp control)*): each recording consists of one or more stimulus blocks that contain onset/offset times and stimulus data that are generated using *RunExperiment.m* and its dependency functions. The stimulation scripts automatically record the time stamps of the national instruments (NI) I/O card and internal stimulus times (PC times) generated by Psychtoolbox. Second column (*SpikeGLX*): although the NI I/O streams and IMEC neuropixels (Npx) data may appear to be monolithic, they are entirely separate data streams with independent clock times. SpikeGLX produces data streams at three different sampling frequencies: the NI I/O data, the AP channels, and the downsampled LFP channels. To synchronize these data streams, you need common pulses; i.e., connect the IMEC pulse generator to an NI I/O stream. Third column (*RunEyeTracker*): optionally, you can use the eye-tracker program to perform online eye-tracking and automatically synchronize your pupil tracking movie to the rest of your data. While recording, the eye tracker program periodically queries SpikeGLX for NI data stream time stamps and saves a log file containing synchronous video time stamps and NI time stamps. Moreover, you can define an ROI for luminance-based high-resolution synchronization (LED pulses). Last column (*Histology*): by registering your histology slices to the Allen Brain Atlas using *AP_histology*, *SHARP-track*, or a similar pipeline, you can extract your probe's approximate location. Using the coordinate adjustment GUI you can then fine-tune the coordinates based on electrophysiological markers, such as responsiveness to experimental treatment and the across-channel spiking correlation. The *RecordingProcessor* allows you to easily cluster your spikes with Kilosort (1), combine and synchronize data from multiple sources (2), adjust probe coordinates (3) and automatically extract which brain area each cell was located, using the ABA API (4).

Installation instructions

Before you start:

1. Make sure your computer is up-to-date, and has compatible (not necessarily the newest) CUDA drivers for your Nvidia GPU. If you don't have an Nvidia GPU, you'll have to buy one, because multiple Acquipix functions, as well as the external libraries it uses (Kilosort, Psychtoolbox) make use of CUDA-specific GPU-accelerated computing.
2. You will want a fast SSD for data buffering for various pre-processing operations. You can run everything without an SSD, but it will be noticeably slower.

Installing the external libraries for the recording computer:

1. If you want to run visual experiments on this computer, download Psychtoolbox and follow the installation instructions: <http://psychtoolbox.org/download>
2. If you want to record Neuropixels on this computer, install SpikeGLX: <https://github.com/billkarsh/SpikeGLX>
3. If you want to perform pupil-tracking specifically, or generally want to have accurate & automatic synchronization between your recorded videos, stimulation and Neuropixels data, you can install this Acquipix-independent module on the computer where you're recording the video's: <https://github.com/JorritMontijn/EyeTracker>

Installing the external libraries for preprocessing your data:

1. Download the Kilosort (v3 or v2.5 are both supported), Npy-matlab and spikes repositories and follow their installation instructions: <https://github.com/MouseLand/Kilosort/>, <https://github.com/kwikteam/npymatlab> and <https://github.com/cortex-lab/spikes>. To compile Kilosort's GPU functions, you'll first need to install a compiler. You can for example use visual studio, but note that only specific versions will work in combination with specific versions of matlab: <https://visualstudio.microsoft.com/vs/older-downloads/>. In all cases, make sure you get the **Community** version, and make sure you also install the **C++ compiler**.
2. Compiling Kilosort's CUDA functions with **Visual Studio Community 2015** with Update 3 seems to work with Matlab R2019b, but the version compatibility seems somewhat arcane and random. If your matlab is R2019b or later, try VSC2015 first. If it doesn't work, you can try a different version (i.e., VSC2013).
3. Rename and edit the config and chanmap files to match your preferred settings and copy them to a folder outside the git repositories if you want to make sure they don't get overwritten accidentally.

Installing Acquipix and its dependencies:

1. Download the Acquipix repository at <https://github.com/JorritMontijn/Acquipix>
2. Download the GeneralAnalysis and MNCP repositories from <https://github.com/JorritMontijn/>
3. Optional: If you wish to automatically compute a responsiveness metric for your putative single cells, download <https://github.com/JorritMontijn/ZETA>

4. Optional: If you wish to save higher quality figures, download <https://github.com/JorritMontijn/exportFig>

Workflow summary:

To run an experiment:

1. Start SpikeGLX and begin a new acquisition
2. Start "RunExperiment" in MATLAB
3. Select your stimulation and parameters & start the experiment (spikeglx will start recording automatically)

Optional steps for pre-processing:

- A) Pre-process your eye-tracking data
- B) Pre-process your probe coordinates

To pre-process your spiking data with six mouse clicks:

1. Start "RunRecordingProcessor" in MATLAB and compile your data library
2. Select the recording you wish to preprocess
3. Cluster your data by clicking the button in the GUI
4. Combine data from multiple sources by clicking another button using the GUI
5. Adjust the probe coordinates and make some changes if necessary
6. Then export with click #6

User guide

Running an experiment

To run an experiment, run matlab and execute “RunExperiment”. This will ask you for the name of your experiment and/or subject. It will automatically add the date to the name when saving the output files. After entering the name, it will show you this screen:

Name: MrMouse

SpikeGLX Client

Host: 192.87.10.238 Idle

Recording: ...

Free space: NaN GB Start Rec. ☒ SGL

NI-DAQ Client

Device:

Pupil LED gain: 1

Sync LED gain: 0.5

LEDs On LEDs Off ☒ Daq

Select Experiment:

RunReceptiveFieldMapping

Select Stimulus Set:

Set 1

strExpType

RunReceptiveFieldMapping

dblSecsBlankAtStart

3

dblSecsBlankPre

0.3

dblSecsStimDur

0.6

dblSecsBlankPost

0.1

dblSecsBlankAtEnd

3

Evaluate Inputs

Estimated Duration:
01:13 (mm:ss); 68x1

Save New Set

Overwrite

Delete

strStimType

'SparseCheckers'

strLinLoc

'LinLoc matrix is upside down; top of screen are bottom elements'

strOutputPath

'C:_Data\Exp'

strTempObjectPath

'X:\WorritMontijn\'

dblSubjectPosX_cm

0

Start Experiment!

The screen consists of three panels, from top to bottom:

- The experiment-independent settings, including the name of the recording (“MrMouse”), the SpikeGLX connection settings (in the above window it’s not connected yet), and the NI-DAQ settings, where you can set two LED Daq outputs. Below you can select the experiment (“RunReceptiveFieldMapping”) and which stimulus preset you wish to use.
- Selecting an experiment and stimulus preset will load the preset settings into the middle panel. You can edit the values here and create a new preset.
- The bottom panel show default parameter values that the program extracted from the experiment’s code; all fields of the structure “sStimParamsSettings” are listed here. You can overwrite them here, but they will not be saved. The idea of this separation is to have easy and uncluttered access to variables you would often change in the stimulus presets (middle panel), while still being able to change all other, less often changed parameters, in the bottom panel.

When you are happy with your selection and wish to run the experiment, you can click “Start Experiment!”. If SpikeGLX and the EyeTracker are not yet recording, they will now start automatically.

Using the online analysis tools

Acquipix comes with a set of online analysis tools that allow you to stream SpikeGLX data in real-time and calculate MUA tuning curves per channel: runOnlineOT for Orientation Tuning, runOnlineNM for Natural Movies, and runOnlineRF for Receptive Field mapping.

Preprocessing your data

Setting the variables and parameters

To start preprocessing, run matlab and execute “runRecordingProcessor”. The startup screen will look like this:

The screenshot shows a MATLAB GUI titled "Master paths". It contains several input fields for setting paths and parameters. The fields are arranged in a list, each with a label and a text box. The labels are: "Set Output:", "Set Temp:", "Set Ephys:", "Set StimLog:", "Set EyeTracking:", "Set ProbeLoc:", "Set AllenCCF:", and "Set ConfigFile:". The corresponding text boxes contain the following paths: "H:\DataPreProcessed", "E:_TempData", "H:\DataNeuropixels", "H:\DataNeuropixels", "H:\DataNeuropixels", "D:\Data\Raw\Histology", "F:\Data\AllenCCF", and "F:\Code\Acquisition\Acquix\subfunctionsPP\configFile384_Npx3B2.m". Below these fields are two buttons: "Edit variables" and "Find data". To the right of these buttons are two labels: "Temp space:" and "Spike sorter:". The "Temp space:" label is followed by the text "146.9 GB". The "Spike sorter:" label is followed by the text "Kilosort2". To the right of these labels is a checkbox labeled "Keep temp". At the bottom of the GUI are five buttons: "Cluster spikes", "Combine sources", "Adjust coords", "Export files", and "Transform to NWB".

To start pre-processing your data, you will first need to set the correct paths:

“Set Output:” specifies where the aggregated and pre-processed data files will be stored.

“Set Temp:” path to fast SSD for data buffering (e.g., used by KiloSort).

“Set Ephys:” root path of your SpikeGLX data files. The actual data files can be in subdirectories. The recommended path structure is something like /root/experiment-set/recordingdate/. To avoid having to compile huge libraries on startup, it is also recommended to use a different root directory for each user and project, but that is of course all up to you.

“Set StimLog:” Same as above, but for the stimulus logs produced by RunExperiment’s subsidiary stimulation functions (e.g., RunDriftingGrating). If you wish, you can put these files together with the SpikeGLX files in the same folder and set the same root. The library compiler combines files from the same recording based on the SpikeGLX run name if they are stored as a variable in the file. If not, it defaults to using the experiment ID in the file name. It is therefore recommended to not change the file names yourself, as this might potentially cause incorrect file attributions.

“Set EyeTracking:” Same as above, but for the output of the EyeTracker repository. Use of the EyeTracker is optional.

“Set ProbeLoc:” Default path to select a probe coordinate file. As there is no way to automatically match the output of AP_histology, SHARP-track, or others to a specific recording, you will have to manually link them by selecting the correct probe coordinate file.

“Set AllenCCF:” Path to the Allen Brain Atlas Common Coordinate Framework data files.

“Set ConfigFile:” Path to the preprocessing configuration file used by Kilosort. A default file for the Neuropixels 1.0 probe is supplied in the /Acquipix/subfunctionsPP/ folder (configFile384_Npx3B2.m). You can edit the file yourself to match your probe properties, but it is recommended to save this file to a separate directory so you don’t make changes to files the /Acquipix/ directory. This way you can easily update the code with github if a new version is released. If you start locally overwriting the repository files, you might have to create a fork and merge commits manually.

Next, you will want to edit the variables using the “Edit variables:” button. This will make another screen appear, where you can set some metadata that you want to be exported as a json file along with your aggregate data file:



The screenshot shows a window titled "Edit variables" with a "Help" button in the top right corner. At the top, there are "Load" and "Save as" buttons, followed by a text field containing the path "D:\Data\Raw\Neuropixels\Acquipix\Metavars\Wi". Below this is a table of variables, each with a "Rem" button to its left. The variables and their values are:

Variable	Value
version	1.0
dataset	Neuropixels_data
investigator	Jorrit_Montijn
project	MontijnNPX2020
setup	Neuropixels
stimulus	VisStimAcquipix
condition	none
subjecttype	DBA
niCh0	onset
niCh1	@PP_GetRunSpeed
niCh2	sync

Below the table is a "New field" button. At the bottom of the window is an "Accept" button.

You can add and remove variables that will be added to the json file as you wish. There is one “special” type of variable here, and that is anything starting with “niCh”. In the above window, niCh0 with the value “onset” specifies that channel 0 is the channel in the National Instruments data file (nidq) that can be used to synchronize stimulus onsets. E.g.: on my setup, we use a photodiode stuck to the upper right corner of the screen where a white rectangle is presented at the beginning of each stimulus presentation. This photodiode then sends a voltage signal to NI channel 0 (index 1 in matlab) so we know exactly when the visual stimulus appeared on screen. It is not required to use

this additional synchronization, but if you don't, there will likely be multiple milliseconds of uncorrected lag+jitter between when your stimulus PC logs the stimulus onset and when it actually appears on the screen.

Likewise, niCh2 has the value "sync", indicating this is the pulse channel receiving 1 second pulses from the IMEC card. You can specify any number of additional National Instruments channels. In the above window, niCh1 is set to "@PP_GetRunSpeed": on our setup, NI channel 1 receives encoder input from a running wheel, which during the synthesis combination step is fed into a data processor function with the name PP_GetRunSpeed, as described in the following section.

Custom data-preprocessors of NI channels

Using the above example, we're feeding the raw binary data collected on NI channel 1 (matlab index 2) into a data pre-processor function PP_GetRunSpeed that transforms voltages into running speeds. This function is called when running the data synthesis combination step. The syntax is the following:

```
output = DataProcessorFunction(vecBinaryChannelData, sMetaNI)
```

The first input (`vecBinaryChannelData`) is the raw NI channel data and the second input is a structure containing the SpikeGLX metadata of the NI file, including for example the NI sampling frequency. The function can produce any output, including structures with subfields, but it will only save the first output into the exported data file. The output will be listed as a field in the `sAP` structure and have the same name as the processor function. In the case of `PP_GetRunSpeed`, this produces a field `sAP.PP_GetRunSpeed`. This field itself is also structure, with the following fields:

`sAP.PP_GetRunSpeed.vecOutT` is a vector with subsampled timestamps (`t0=0`) (in seconds)

`sAP.PP_GetRunSpeed.vecTraversed_m` is the traversed distance per time step (in meters)

`sAP.PP_GetRunSpeed.vecSpeed_mps` is the filtered running speed over the last second (in m/s)

Compiling the library and pre-processing the data

Now that all parameters have been set, you can compile the data library by clicking "Compile data library". The program will now compile a list of all your data files and group them by experiment. The leading file is the SpikeGLX nidq meta file:

Master paths

Set Output: H:\DataPreProcessed

Set Temp: E:_TempData

Set Ephys: H:\DataNeuropixels

Set StimLog: H:\DataNeuropixels

Set EyeTracking: H:\DataNeuropixels

Set ProbeLoc: D:\Data\Raw\Histology

Set AllenCCF: F:\Data\AllenCCF

Set ConfigFile: F:\Code\Acquisition\Acquix\subfunctionsPP\configFile384_Npx3B2.m

Edit variables Find data Temp space: 146.9 GB Spike sorter: Kilosort2 ☐ Keep temp

Run?	Synthesized	Clustered	Pupil data	Probe coords	Slim files	AP Ephys	LFP Ephys	NI Ephys	File name	Date	Total size
<input type="checkbox"/>	Y	Y	1	3	4	Y	Y	Y	RecMA5_2021-02-25R01_g0_t0	20210225	128.8GB
<input type="checkbox"/>	N	Y	1	1	3	Y	Y	Y	RecMA52_2021-02-26R01_g0_t0	20210226	66.5GB
<input type="checkbox"/>	Y	Y	1	1	3	Y	Y	Y	RecMA5_2021-02-26R01_g0_t0	20210226	80.3GB
<input type="checkbox"/>	Y	Y	1	4	4	Y	Y	Y	RecMA5_2021-03-01R01_g0_t0	20210301	107.4GB
<input type="checkbox"/>	Y	Y	0	1	0	Y	Y	Y	Runtest_2021-06-11_g0_t0	20210611	1.1GB

Cluster spikes Combine sources Adjust coords Export files Transform to NWB

This will show you a list of all experiments, their dates and sizes, and whether certain (preprocessed) files are present or not. Hovering over the different indicators will give you further detailed tooltip information.

You can now select multiple recordings (mark the check box under the "Run?" column) and click "Cluster spikes", after which the program will run the indicated spike sorter that it found on your matlab path on all your selected files in serial. When the clustering is finished, you can combine all files of one recording and synchronize the timestamps of the various data into a single synthesized aggregate by clicking "Combine sources". Now you can link the probe coordinate file by clicking on the button under "Probe coords", which will give you another pop-up showing the probe number and area of entry of all probes listed in this file. Select the correct probe and click OK:

Select probe # for RecMA5_2021-03-01R01_t0

Probe 1, starting at "Primary visual area layer 1"

Probe 2, starting at "posteromedial visual area layer 1"

Probe 3, starting at "Anteromedial visual area layer 1"

OK Cancel

Now that all data has been combined, you can fine-tune the probe location by clicking “Adjust coords”. This will open a new window (which you can also access by running `getProbeHistology.m`):



- 1) Similar to the SHARP-track brain browser, it shows you an outline of the mouse brain, the current location of your probe (blue line), the probe's entry point into the brain (red cross), the histological locations you selected for this probe (blue dots), and a transparent slice through the brain showing the location of the ABA-annotated brain areas.
- 2) The plot on the upper right shows the responsiveness of your clusters (blue dots) as a function of depth (y-axis) and ZETA score (x-axis). The higher the score, the more likely your cell is responding to any of your stimuli (see <https://doi.org/10.1101/2020.09.24.311118> for more information on the ZETA-test). The grey lines are area boundaries along the probe's trajectory, as shown on the right. Note that based only on ZETA scores the upper and lower boundaries of visual cortex (area PM) are already very clear. The same goes for when the visual nuclei in the midbrain appear, starting with area NOT in the above example.
- 3) The MUA correlation plot shows how correlated the spiking activity is between different channels on the probe. Usually channels that are in the same area will show a higher correlation with each other than with other channels. The MUA plot is especially useful for determining the location of the pia mater, as channels outside the brain are often highly correlated due to shared noise.
- 4) Finally, the template depth & rate plot shows the normalized firing rate of all clusters (i.e., kilosort templates) as a function of depth. Certain areas may show higher firing rates than others, which can give you another clue of where the area boundaries might lie.

When you are done, you can close the GUI and click “save and exit” to save the adjusted coordinates.

Exporting data to the Acquipix (AP) structure format

Once you finished the above steps, you can click “export files” in the RecordingProcessor to assign areas to each cluster, save an aggregate recording file, and export a json file with metadata to your output directory. This creates a file containing an sAP structure. In short, the sAP structure contains a cell array with stimulation information per block (cellBlock), data on each cluster, such as its spike times and brain area (sCluster), the source data (sSources), json data (sJson), probe location, and optionally pupil tracking information and any other data structures generated from additional NI channels. See the section “***The Acquipix (sAP) data structure format***” below for more detailed information on the sAP structure.

Exporting data to the Neurodata Without Borders (NWB) format

If you wish to share your data or simply prefer to work with the NWB format, you can also export the data to the NWB data format by clicking “Transform to NWB”. (to do...)

Technical descriptions

If you want to make your own stimulation and/or online analysis scripts, you can use the modules present in Acquipix. I've tried to make everything as modular as possible, so it should be easy to plug in experiments of your own design.

The Acquipix (sAP) data structure format

Overview

The sAP structure is subdivided into several substructures. Let's look at an example recording:

```
sAP =  
  
    struct with fields:  
  
        cellBlock: {[1×1 struct] [1×1 struct] [1×1 struct]}  
        sPupil: [1×1 struct]  
        sCluster: [1×461 struct]  
        sSources: [1×1 struct]  
        sJson: [1×1 struct]  
        vecProbeCoords: [-3289 1914 2790 11.9771 22.5562]  
        ProbeCoordsDesc: 'AP, ML, depth, AP-angle, ML-angle'
```

The above shows the coordinates of brain entry relative to Bregma, the depth of the tip and AP and ML angles. It contains a cell array `cellBlock`, where each element is a structure containing data on a single recording block, for example a set of drifting grating presentations (see below). It also always contains the structure `sCluster`, where each element corresponds to a single cluster (SU, MU, or noise) that was defined using the spike sorter in the clustering step (see further below). The `sJson` structure contains the same data as the `.json` file; and if you performed eye-tracking you will also have a structure called `sPupil`. In addition, the `sAP` structure may contain miscellaneous structures containing data on custom pre-processed NI data streams if you specified any in the meta variable list.

The cellBlock structure

An example `cellBlock` structure looks like this:

```
sAP.cellBlock{2} =  
  
    struct with fields:  
  
        strExpType: 'RunDriftingGratings'  
        strSyncType: 'Good: NI timestamps'  
        intNumRepeats: 20  
        vecOrientationNoise: 0  
        intTrialNum: 480  
        intStimNumber: 480  
        dblStimFrameDur: 0.0167  
        SampRateNI: 8.4745e+04  
        vecTrialStimTypes: [1×480 double]  
        TrialNumber: [1×480 double]  
        OrientationNoise: [1×480 double]  
        Phase: [1×480 double]  
        vecStimOnTime: [1×480 double]  
        vecStimOffTime: [1×480 double]
```

```

vecPupilStimOnTime: [1×480 double]
vecPupilStimOffTime: [1×480 double]
sStimParams: [1×1 struct]
sStimObject: [1×24 struct]

```

It lists the type of experiment (here: drifting gratings), whether the synchronization could use national instruments timestamps or not and some optional variables such as the stimulus times in pupil video timestamps, onset phase, and the number of trials. However, the most important variables here are the following four:

- `vecStimOnTime`: synchronized stimulus onset times aligned to IMEC spike times
- `vecStimOffTime`: same as above, for the offset
- `vecTrialStimTypes`: the stimulus type with which you can access the stimulus parameters using `sStimObject`
- `sStimObject`: a structure containing the parameters for each unique stimulus type.

For example, to retrieve the stimulus onsets, offsets, and parameters for the sixth trial, you can do this:

```

sBlock = sAP.cellBlock{2};
Trial = 6;
StimOnset = sBlock.vecStimOnTime(Trial);
StimOffset = sBlock.vecStimOffTime(Trial);
StimType = sBlock.vecTrialStimTypes(Trial);
TrialObject = sBlock.sStimObject(StimType);

```

The sStimObject substructure

The stimulus parameters for each unique stimulus can be found in `sStimObject`, and the fields will depend on the type of stimulus. For example, for a drifting grating stimulus we can find:

```

sAP.cellBlock{2}.sStimObject(6) =

struct with fields:

    StimType: 'SquareGrating'
    CornerTrigger: 2
    CornerSize: 0.0333
    ScreenDistance_cm: 17
    SubjectPosX_cm: 0
    SubjectPosY_cm: -2.5000
    AntiAlias: 0
    UseGPU: 1
    StimPosX_deg: 0
    StimPosY_deg: 0
    StimulusSize_deg: 140
    SoftEdge_deg: 2
    Background: 0.5000
    UseMask: 0
    Contrast: 100
    Luminance: 100
    Orientation: 75
    OrientationNoise: 0
    SpatialFrequency: 0.0500
    TemporalFrequency: 1
    Phase: NaN
    FrameRate: 60

```

This shows a whole list of different parameters, but you will presumably be most interested in the orientation in case of drifting gratings (here shown in degrees). Note that the starting phase is randomized stochastically, so this is not saved in the `sStimObject` structure. Rather, if you wish to access the starting phase of the stimulus, you can find it in `sAP.cellBlock{2}.Phase`.

The sCluster structure

Perhaps the most important information of your recording can be found in `sCluster`, which lists information on each cluster in your recording, including the brain where it was found, the spike times, and – if you have the ZETA repository installed – the responsiveness p-value for the different stimulus blocks. The first cluster in the example recording looks like this:

```
sAP.sCluster(1) =  
  
    struct with fields:  
  
                Exp: 'RecMA8_2021-02-20R01_g0_t0'  
                Rec: '20210220'  
                Area: 'Midbrain'  
    SubjectType: 'DBA'  
        Subject: 'RecMA8_  
        Date: '2021-08-02'  
        Depth: 3.8000e+03  
        Cluster: 1  
        IdxClust: 0  
        SpikeTimes: [9115×1 double]  
    NonStationarity: 0.0098  
    Violations1ms: 0.1289  
    Violations2ms: 0.1741  
    Contamination: 100  
    KilosortGood: 1  
                ZetaP: [0.9205 0.9181 0.9874]  
                MeanP: [0.6184 0.5942 0.4306]  
    DepthBelowIntersect: 2750
```

We can see which experiment this cluster is from, the Allen Brain Atlas area name it was found in (Area: 'Midbrain'), the spike times (SpikeTimes), and some information on whether this is a good cluster or not. This includes the following variables:

- **NonStationarity:** this value indicates whether the spikes of this unit are found exclusively in the beginning of the recording (value of -1), exclusively at the end (value of +1), or equally distributed between beginning and end (value of 0). This cell shows a very low non-stationarity of 0.0098. Note, however, that if 50% of the spikes occur in the first second of the recording, and the other 50% during the last second, this value will still be 0. So the interpretation here is: values close to -1 or +1 are bad; values close to 0 mean that the cell *might* have a low temporal non-stationarity across the recording.
- **Violations1ms:** normalized number of spikes occurring within 1ms of another spike. If you would shuffle the spike times randomly and recompute this value, this would give you a value of 1. Values below 1 mean that this unit shows a refractory period. The closer to 0 the better. This is an indication of the spike contamination calculated somewhat differently from the kilosort version (see below).
- **Violations2ms:** same as above for a 2-ms window around a spike.
- **Contamination:** the estimated percentage of spikes that belong to other clusters that are contaminating this unit. Unfortunately the current version of Kilosort (as of August 2021) has

an error where it can erroneously give high contamination values, so take this value with a grain of salt. Especially if it is exactly 100 you can safely ignore this value, since a contamination of exactly 100 is theoretically impossible.

- KilosortGood: this value is 1 if Kilosort determined this was a good cell; 0 otherwise.

Experiment scripts

Overview

We will use “RunDriftingGratings.m” as the example here, but other files work the same. All experiment scripts are built so they can be run as a script in the base workspace. The idea here is that your variables will be available in the base workspace even in the case you get an error, forget to save your data, there’s a network drive malfunction, etc. I’ve also built the experiment script so it can be run both independently with default parameters, and as a subfunction of the RunExperiment GUI with supplied inputs.

Description of inputs and structure of an experiment script

When you press “Start Experiment!” in RunExperiment, it runs the following code block:

```
%run!
assignin('base','sExpMeta',sExpMeta);
assignin('base','sStimParamsSettings',sStimParamsSettings);
assignin('base','sStimPresets',sStimPresets);

evalin('base',strStimType);
```

As mentioned above, you can see that this assigns three variables to the base workspace: sExpMeta, sStimParamsSettings, and sStimPresets. The experiment you wish to run is contained in strStimType; for example runDriftingGratings.

sExpMeta has fixed fields filled by RE_genGUI:

```
sExpMeta.boolUseSGL = boolUseSGL;
sExpMeta.boolUseNI = boolUseNI;
sExpMeta.dblPupilLightMultiplier = dblPupilLightMultiplier;
sExpMeta.dblSyncLightMultiplier = dblSyncLightMultiplier;
if boolUseSGL
    sExpMeta.strHostAddress = sRE.strHostAddress;
    sExpMeta.hSGL = sRE.hSGL;
    sExpMeta.sParamsSGL = sRE.sParamsSGL;
    sExpMeta.strRunName = sRE.strRunName;
end
if boolUseNI
    sExpMeta.objDaqOut = sRE.objDaqOut;
end
```

It contains some required static information about this recording; whether you wish to use SpikeGLX (boolUseSGL), the National Instruments I/O box (boolUseNI), gain values for the pupil and synchronization LEDs. If SpikeGLX is connected, it supplies the server’s IP address (strHostAddress), a handle to the SpikeGLX object (hSGL), a parameter structure supplied by

SpikeGLX (`sParamsSGL`), the name of the current recording (`strRunName`); and if the NI I/O box is used, it supplies a Matlab Daq Interface object `objDaqOut`.

`sStimParamsSettings` is a structure with parameters that are dynamically extracted from the experiment script code itself. If you look in `runDriftingGratings.m` you can see the following block:

```
%% input params
fprintf('Loading settings...\n');
if ~exist('sStimParamsSettings','var') || isempty(sStimParamsSettings)
    %general
    sStimParamsSettings = struct;
    sStimParamsSettings.strStimType = 'SquareGrating';
    sStimParamsSettings.strOutputPath = 'C:\_Data\Exp'; %appends date

    %visual space parameters
    sStimParamsSettings.dblSubjectPosX_cm = 0; %cm; relative to c of scr
    sStimParamsSettings.dblSubjectPosY_cm = -2.5;%cm; rel to c of scr

    (...)
else
    % evaluate and assign pre-defined values to structure
    cellFields = fieldnames(sStimParamsSettings);
    for intField=1:numel(cellFields)
        try
            sStimParamsSettings.(cellFields{intField}) = ...
                eval(sStimParamsSettings.(cellFields{intField}));
        catch
            sStimParamsSettings.(cellFields{intField}) = ...
                sStimParamsSettings.(cellFields{intField});
        end
    end
end
```

The file `RE_getParams.m` reads the experiment script file and loads the field name, value and comment as tooltip information. Any default value you put/change in your experiment file will therefore automatically pop up in the `RunExperiment` GUI. Moreover, if you change values in the GUI, it will supply them in the `sStimParamsSettings` structure and load them in the `else` block. These are usually semi-static variables, like the experiment script's stimulus type (i.e., "SquareGrating"), the location of the mouse relative to the screen, etc. Finally, `sStimParamsSettings.strRecording` is always added by `RunExperiment` to identify the current recording.

`sStimPresets` has dynamic parameter values that you might want to change often and save as a separate pre-set; e.g., the number of repetitions, orientations, stimulus duration, etc. You can change these values easily in the GUI, and on first run, it generates pre-set files in `assertPresets.m`. For example, you can see the following block here:

```
elseif strcmp(strExp, 'RunDriftingGratings')
    for intSet=1:3
        if intSet == 1
            sStimPresets = struct;
            sStimPresets.strExpType = strExp;
            sStimPresets.intNumRepeats = 100;
            sStimPresets.vecOrientations = 0:15:345;
```

```

sStimPresets.vecOrientationNoise = 0;
sStimPresets.dblSecsBlankAtStart = 3;
sStimPresets.dblSecsBlankPre = 0.4000;
sStimPresets.dblSecsStimDur = 1;
sStimPresets.dblSecsBlankPost = 0.1000;
sStimPresets.dblSecsBlankAtEnd = 3;
(...)

```

If you wish to make your own experiment file, you will have to edit `assertPresets` and add your own default values here. Note that you can make it into anything you wish here, with the exception of the `strExpType` field: this is compulsory to differentiate the different experiment scripts. It must match the name of your experiment. For example, because the experiment script is called `RunDriftingGratings.m`, `strExpType` takes the value `'RunDriftingGratings'`.

Structure of the experiment script & output file

You can basically do anything here, as long as you observe the restrictions imposed by the three input variables described above, and the structure of the output file. Looking in `RunDriftingGratings`, we can see that the output file is saved as:

```
save(fullfile(strLogDir, strFilename), 'structEP', 'sParamsSGL');
```

The `sParamsSGL` structure is the same as supplied by `sExpMeta`, so you can use the following:

```
sParamsSGL = sExpMeta.sParamsSGL;
```

The `structEP` structure has four obligatory fields:

`structEP.ActOnSecs` is the actual stimulus onset time in seconds (e.g., the PsychToolBox Screen flip time of the first frame of your stimulus)

`structEP.ActOffSecs` is the actual stimulus offset time in seconds (e.g., the PsychToolBox Screen flip time of the first blanking frame after your stimulus ended)

`structEP.ActOnNI` is the actual stimulus onset in National Instruments (nidq) time stamp

`structEP.ActOffNI` is the actual stimulus offset in National Instruments (nidq) time stamp

You can obtain the National Instruments time stamp using:

```

hSGL = sExpMeta.hSGL;
intStreamNI = -1;
dblSampFreqNI = GetSampleRate(hSGL, intStreamNI);
dblTimeStampNI = GetScanCount(hSGL, intStreamNI)/dblSampFreqNI;

```

That said, it is of course wise to save more information than that; in the case of drifting gratings this would be the orientation during each trial, the spatial frequency, etc. To ensure I always have all meta information, I also include the `sStimParams` structure before saving:

```

structEP.sStimParams = sStimParams;
structEP.sStimObject = sStimObject;
structEP.sStimTypeList = sStimTypeList;
save(fullfile(strLogDir, strFilename), 'structEP', 'sParamsSGL');

```

But of course you can decide differently for yourself. A minimal experiment script could therefore look something like this:

```
% evaluate and assign pre-defined values to structure
if ~strcmpi(sStimParamsSettings.strStimType, 'StupidExperiment')
    sStimParamsSettings = struct;
    sStimParamsSettings.strStimType = 'StupidExperiment';
    sStimParamsSettings.strOutputPath = 'C:\Data\';
    sStimParamsSettings.strTempObjectPath = 'C:\Temp\';
else
    cellFields = fieldnames(sStimParamsSettings);
    for intField=1:numel(cellFields)
        try
            sStimParamsSettings.(cellFields{intField}) = ...
                eval(sStimParamsSettings.(cellFields{intField}));
        catch
            sStimParamsSettings.(cellFields{intField}) = ...
                sStimParamsSettings.(cellFields{intField});
        end
    end
end

% set output locations for logs
strRecording = sStimParamsSettings.strRecording;
strOutputPath = sStimParamsSettings.strOutputPath;
strTempPath = sStimParamsSettings.strTempObjectPath;
strThisFilePath = mfilename('fullpath');
[strFilename, strLogDir, strTempDir, strTexDir] = ...
    RE_assertPaths(strOutputPath, strRecording, strTempPath, strThisFilePath);

%get SGL data
hSGL = sExpMeta.hSGL;
strRunName = sExpMeta.strRunName;
sParamsSGL = sExpMeta.sParamsSGL;
intStreamNI = -1;
dblSampFreqNI = GetSampleRate(hSGL, intStreamNI);

%run trials
structEP = struct;
hStartTic = tic;
for intThisTrial = 1:sStimPresets.intNumRepeats
    %"inter-trial" interval
    pause(0.5);
    structEP.ActOnSecs(intThisTrial) = toc(hStartTic);
    structEP.ActOnNI(intThisTrial) = GetScanCount(hSGL, ...
        intStreamNI)/dblSampFreqNI;

    %"stimulus" period
    pause(1);
    structEP.ActOffSecs(intThisTrial) = toc(hStartTic);
    structEP.ActOffNI(intThisTrial) = GetScanCount(hSGL, ...
        intStreamNI)/dblSampFreqNI;
end

%save data
save(fullfile(strLogDir, strFilename), 'structEP', 'sParamsSGL');
```

Integrating your experiment file into the Acquipix RunExperiment GUI

If you've followed the descriptions above, the only thing you need to do now is edit the file `assertPresets.m` and add the default parameters for a stimulus preset.

Optional: you can add your experiment to `RE_evaluateStimPresets.m`, which calculates how long your experiment will take based on your preset parameters and ensures there are no incompatible values. It will still work if you don't, but the GUI will show "Cannot resolve duration" when calculating the expected duration of your experiment.

Online analysis functions

Overview

Acquipix comes with three real-time data stream analysis GUIs: `runOnlineOT` for orientation tuning with drifting gratings, `runOnlineNM` for natural movie response analysis, and `runOnlineRF` for receptive field mapping. We will use "runOnlineOT" as the example.

To be continued...

Troubleshooting

Question (“actually, it’s more of a comment”): *It doesn’t work*

Answer: Restart your PC

Q: *I downloaded everything, but it says files are missing*

A: Double check you have added all folders to the path in Matlab, you have the required Matlab toolboxes installed (Curve Fitting, Parallel Computing), and you’re using a supported matlab version: R2019b is tested and works; anything earlier than R2016b will fail for sure; other versions might work. You may also need additional toolboxes, such as Image Processing and Acquisition to perform eye-tracking. If it’s still not working after you’ve tried the above, google the filename and reinstall its source repository. If it still fails, create a report here:

<https://github.com/JorritMontijn/Acquipix/issues>.

Q: *I cannot compile Kilosort’s GPU code*

A: First try installing VSC2015 (make sure you have the Visual Studio **Community** version) and make sure you install the C++ compiler. Then **restart your PC** and try again. If it doesn’t work: uninstall and try VSC2013. If neither of them work, look for help here: <https://github.com/MouseLand/Kilosort/>

Q: *I cannot run any GPU code in matlab (i.e., gpuArray() fails)*

A: Make sure that you have the correct CUDA drivers installed for your GPU. Note that if you’re using anything other than a (modern) Nvidia GPU, you cannot run CUDA.

Q: *I found a bug*

A: Great! Or at least, it’s great that you found it, not that it’s there. If you’ve fixed it, you can make a pull request, otherwise you can create a bug report here:

<https://github.com/JorritMontijn/Acquipix/issues>. Please copy/paste the matlab error message and as much detail as you can about what you were doing when it happened. If I cannot recreate the issue, I probably won’t be able to fix it.